

Uplift Commode Assist:

The lifting commode chair that gives independence and privacy when it's needed most.

Uplift Commode Assist is a self-powered, lifting commode chair that helps users gently seat and raise themselves. It provides safe, controlled support and lifts up to 70% of a user's weight.

Versatile and sturdy, the Commode Assist can be used standalone or over the toilet. For safety, it automatically locks into position when fully depressed.

Benefits

- Helps keep people active and independent.
- Relieves stress and strain on joints.
- Lifts only as needed, allowing users to maintain muscular strength.

Who can benefit

- People who have lost sufficient upper or lower body strength to get up unassisted.
- Those who have lost mobility as a result of:
 - Muscular dystrophy
 - Lou Gehrig's disease (ALS)
 - Post-operative conditions
 - Degenerative joint disease
 - Arthritis
 - Parkinson's disease



Bringing quality to life.






The Commode Assist eases you gently into a seated position and automatically locks in place. To stand, simply press down on the seat release lever (see below).



Using both hands to support yourself, simply shift forward and begin to stand.



The lifting mechanism activates automatically and provides a gentle and stable lift up from the commode.

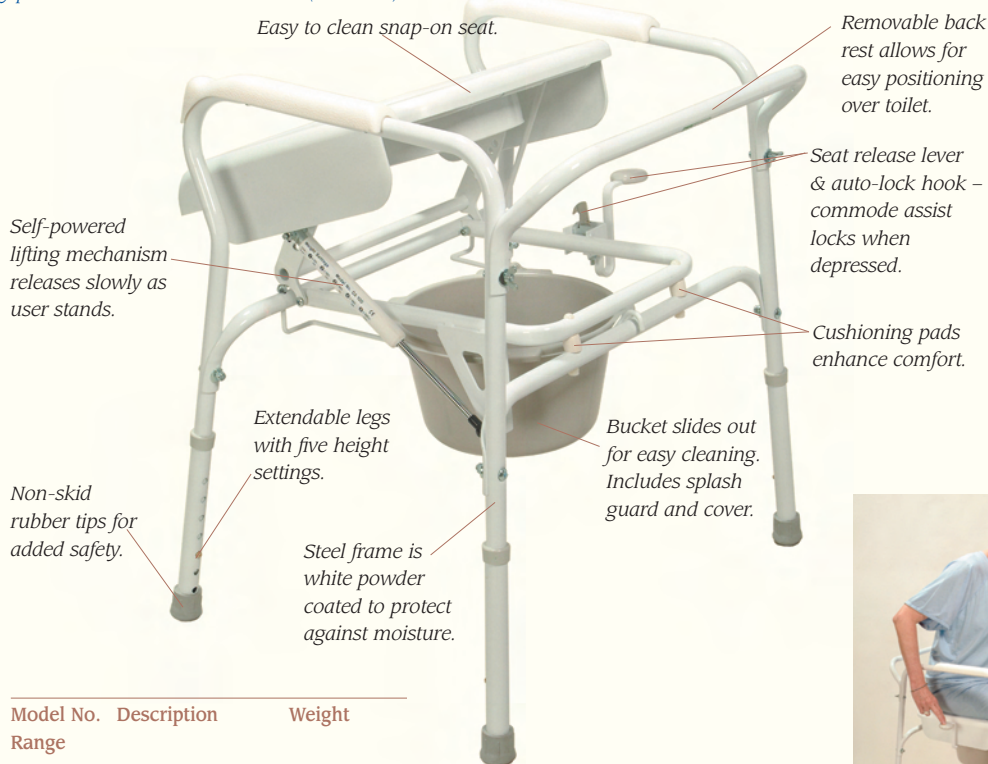

Commode Assist Specifications

Base Dimensions	Depth:	21"/53 cm
	Width:	26"/66 cm
Seat Dimensions	Depth:	16.75"/43 cm
	Width:	16.5"/42 cm
Seat Height	Minimum:	17.75"/45 cm
	Maximum:	23.5"/60 cm
Lift		9"/23 cm @ 20°
Weight		21 lb/9.5 kg



To release the auto-lock, simply push down on the seat release lever. Now both hands are free to use for support when rising (see above sequence).



Model No.	Description	Weight
CA 200*	Lifting commode chair	80-300+ lb 35-136+ kg

*some assembly required



Bringing quality to life.