


The Safety Within

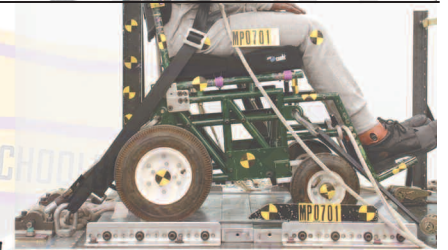

consumer safety QUICKIE® TRANSIT OPTION

QUICKIE®
WHEELCHAIRS

Sunrise Medical is dedicated to making our wheelchairs safe and functional for the consumer. In order to help safeguard individuals who use a wheelchair as their primary seating in a motor vehicle, all Sunrise Medical Transit Option seating is tested to ANSI/RESNA and ISO testing standards.

Quickie Transit Option



What is the Transit Option?

- The Transit Option provides four distinct designated tie-down securement locations on the frame of the wheelchair. Transit Option seating may also include tested mounting hardware used to secure the seating system to the frame of the wheelchair.
- Crash-tested for use with Wheelchair Tiedown and Restraint Systems (WTORS) that comply with SAE J2249. When used in conjunction with proper tie-down and occupant restraints, in appropriate vehicle solutions, the Transit Option offers a complete and approved transportation method.
- Factory installation and Retrofit Kits are available.

What is SAE J2249?

- SAE (Society of Automotive Engineers) created the standards for the WTORS systems and requires that these systems withstand 30 mph crash tests.
- http://www.ercwts.org/RERC_WTS2_KT/RERC_WTS2_KT_Stand/SAE_Restraints_RefDocs/J2249Guide_1_Intro.pdf

What is WC-19?

- WC-19 is a voluntary national US standard called Section 19 ANSI/RESNA WE Volume 1. Sunrise Medical has adhered with this voluntary standard prior to its adoption by the industry in 2000.
- <http://www.ercwts.org/WC19.html>



QUICKIE®
WHEELCHAIRS

WC-19 Tested Products from Sunrise Medical

Pediatric Manual Mobility				Pediatric Power Mobility		Adult Manual Mobility			Adult Power Mobility	

To learn more, contact your local dealer,
call Sunrise Medical at **1-800-333-4000**
or visit **www.sunrisemedical.com**.

* Transit tested in both front and rear wheel drive configurations

* Quickie Groove

© 2008 Sunrise Medical Inc. Sunrise Medical and Quickie are registered trademarks of Sunrise Medical.

930439 Rev. F