


zippie®

GS™


*The ideal system for
the growing pediatric user*

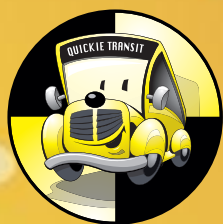

The Zippie GS, named for its built-in Growth System, meets the client's needs in both folding and rigid frames as well as standard and reverse wheel configurations.

Ability to grow the frame as the child's needs change.

The Built-in Growth System easily "grows". The dual seat rails allows for the seating system and accessories to both be mounted without interference in the same small space.

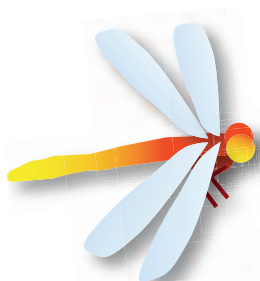

A fully adjustable, universal axle plate allows the GS to be set up in standard or **reverse wheel configurations** and provides half-inch vertical and horizontal adjustments for increased wheel access, without sacrificing growth potential.

Optional One-Arm Drive with Dual Handrims is a proven design that's easy to maneuver and offers trouble-free durability.



ZIPPIE GS SPECIFICATIONS

	Folding	Rigid
HCPCS Code:	E1238	E1237
Overall Weight:*	25 lbs. w/o footrest	29 lbs. w/o footrest
Weight Capacity:	165 lbs.	165 lbs.
Overall Product Width:	19" -27"	19" -27"
Transit Option:	Yes	Yes
Front Seat Height:	14" - 20.25"	14" - 20.25"
Rear Seat Height:	12.5" - 20"	12.5" - 20"
Hanger Options:	60°, 70°, 80°, 90°(Swing-in/out & HD), ELR, ALR	60°, 70°, 80°, 90°(Swing-in/out & HD), ELR, ALR
True Fit Growth System:	Yes	Yes
Warranty:	Lifetime on frame & cross brace, 1 year on parts	Lifetime on frame & cross brace, 1 year on parts

* Actual measurements may vary according to the wheelchair's configuration and accessories