



THANK YOU FOR CHOOSING LUMEX!

Since 1946, Lumex® has been producing quality healthcare products for use in hospitals, long term care facilities, and the home. Lumex has been the innovator in Specialty Healthcare Seating. Often imitated, the Lumex line of Specialty Healthcare Seating continues to set the standard for durability, safety, patient comfort, and ease of use by caregivers.

All of our Specialty Healthcare Seating is backed by our Limited Wraparound Warranty:

- ☐ Frame: Lifetime Warranty
- ☐ Mechanical Components: Three Year Warranty
- ☐ Upholstered Components: One Year Warranty

Please see the GF Health Products, Inc. Warranty for complete details and limitations.

Lumex Specialty Healthcare Seating is designed for use in a variety of environments:

- □ Dialysis
- ☐ Oncology/Drug Delivery
- Ophthalmology
- □ OB/GYN Clinics
- ☐ Same Day Surgery
- ☐ Post-Operative Care
- Bariatric
- Outpatient Care
- □ Blood Collection
- □ Extended Care

Our wide selection of Healthcare Seating can help liberate almost anyone from the confines of bed and make them more comfortable for longer periods of time.

Our line of Clinical Recliners makes care, treatment, and recovery easier and more comfortable. We make our seating modifiable with options to meet special needs. We also make it easy to clean, easy to maintain, and easy to move. We offer recliners, chairs, and more in a variety of attractive styles, colors, and fabrics. And our range of price points ensures that these products fit as neatly into your budget as they do into your facility.

For the most up-to-date information on all G-F products, go to **www.grahamfield.com.**

LUMEX RECLINERS NOW DESIGNED, ENGINEERED, AND ASSEMBLED IN THE U.S.A.

Lumex Specialty Healthcare Seating is manufactured and assembled in our new state of the art facility in Jacksonville, Florida. This new facility allows Lumex the ability to offer you greater flexibility in color and fabric selection, improved quality combined with additional features and options.



LUMEX CLINICAL CARE RECLINERS

All Lumex Clinical Care Recliners have built-in durability and versatility to meet the demanding needs of clinical and critical care use. Lumex Clinical Care Recliners have been designed with input from both patients and caregivers to provide the ultimate in safety, comfort and convenience.



48 COLORS AND PATTERNS

Our FR566, FR577 and FR587 Clinical Care Recliners are now available in a palette of 48 colors! As Interior Designers and Facility Managers play an increasing role in the selection of healthcare seating, the look, fit and finish of the seating used in facilities becomes an integral part of the selection process. Lumex now offers a wide variety of upholstery options and colors. From Sisal to Perriwinkle, Lumex Clinical Care Recliners will fit into any décor or designer specification! Lumex Clinical Care Recliners are available in:

Patriot Plus: 20 Colors in the "Upholstery Standard" of the Healthcare Industry.

Savoy: 10 Colors in a slightly softer hand upholstery.

Renaissanc: 8 Colors in a beautifully patterned upholstery.


Silvertex: **10 Colors** that feature SILVER**GUARD**™ Silver Ion Technology.

Trim thread and welting all match the upholstery color for a smooth, polished look.

PERMABLOK^{3®} PROTECTION!

All Lumex Specialty Healthcare Seating Upholstery features PERMABLOK³ Advanced Vinyl Protective Coating. Engineered to create a tough effective barrier against the three biggest problems encountered in healthcare environments – germs, abrasion, and stains. You can't see it or feel it, yet the protectionit offers results in extended product performance and lasting beauty with a minimum of care. It guards against the surface growth of fungus, mold and mildew spores, which can cause unpleasant odors, unsightly pink and black stains, and even allergic reactions.

All Lumex Upholstery Selections not only keep their good looks but also remain highly resistant to deterioration.

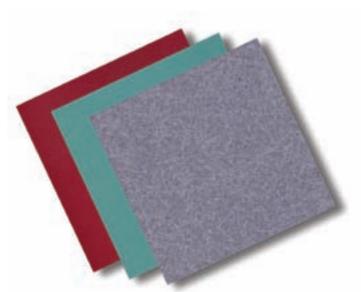

SILVERGUARD® SILVER ION TECHNOLOGY!

Silvertex with SILVER**GUARD** Silver Ion Technology provides a natural, environmentally-friendly, germ-fighting defense to help safeguard against odor and stain-causing bacteria. Silvertex SILVER**GUARD** Silver Ion Technology is EPA registered and FDA approved. When combined with PERMABLOK3, it provides the ultimate performance solution.

How SILVERGUARD® WORKS!

Silver Ions destroy microorganisms by attaching to the bacteria's cell wall to disrupt protein metabolism within the cell. These are slowly released from their inorganic matrix in the presence of moisture. Simply stated—wipe down the upholstery with a damp cloth to activate the Silver Ions!

Silver encapsulates MRSA bacterium and it dies. SILVER**GUARD** is released as needed to effectively control bacterial growth, while leaving sufficient quantities of SILVER**GUARD** to retain antimicrobial protection.



48 Colors Now Available on FR566G, FR566DG, FR577RG & FR587W Clinical Care Recliners

Patriot Plus



Savoy



Renaissance



Silvertex



SILVER**GUARD**®

ALL Lumex Clinical Care Recliners Meet CA TB133!

Safety is of primary importance in all healthcare facilities. All Lumex Clinical Care Recliners meet California Technical Bulletin 133 flammability standards. By utilizing a special flame-blocker material laminated directly onto the upholstery material, Lumex Clinical Care Recliners are not just made from "fire rated" material, but withstand the toughest fire test in use today for healthcare products.

The flexible flame-blocker material does not interfere with the comfort of the recliners because it is laminated directly to the upholstery.

The New Laminated Side Panels not only enhance the look and comfort of the recliner, but feature the same flame-blocker material!

LIQUICELL PRESSURE REDUCTION TECHNOLOGY!

Lumex Clinical Care Recliners are equipped with LiquiCell. LiquiCell is an ultra-thin, liquid-filled interface specifically designed and proven to significantly reduce shear stress, aiding in the prevention of pressure sores and general discomfort.







IMPROVED ARMREST PROFILE

Featuring specially-contoured high-density foam, the new armrest profile enhances patient comfort and support.



SMOOTH SIDE TABLES!

To make cleaning and disinfecting easier, Lumex Clinical Care Recliners feature completely smooth side tables with no seams. The side tables also feature a rolled edge to ensure that prep materials stay in place and not trap contaminants. Side tables feature a simple-to-use mechanism that clearly indicates when table is in the locked position.







Shown in down position

Shown in up position

ULTRAFLEX SUSPENSION SYSTEM—BUILT IN COMFORT!

Patient comfort is a key concern, particularly during long treatment sessions. Lumex Clinical Care Recliners utilize a unique suspension system featuring Ultraflex Tension Bands.

The FR566 Series of Clinical Care Recliners features a series of Ultraflex Tension Bands in the Lumbar Region that support specially contoured foam to ensure proper support and comfort.



OPTIONAL HEAT AND MASSAGE


Our FR566, FR577 and FR587 Clinical Care Recliners are now available with Optional Heat and Massage.

You want your patients to be as comfortable as possible before, during, and after treatment. One of the most common complaints from patients is that they feel cold during treatment. Lumex Clinical Care Recliners are now available with an optional Heat and Massage feature to help your patients stay warm and comfortable! Employing an easy-to-use single switch to independently control both the Heat and Massage feature, the switch is conveniently mounted on the right side of the recliner for easy access / use by both patient and caregiver.

Two heating elements, one in the seat and one in the back, provide gentle heat for the patient. Three Heat Settings allow patients to control just how much warmth they need. The heat feature has an automatic shut-off to help ensure that the patient is always in control of their comfort. Heat automatically shuts off after one hour of continuous use at the same level.

Many patients appreciate a gentle massage during treatment. The massage function offers a four massage settings that control the intensity and location of the massage to ensure that the patient receives the most appropriate level of massage. Massage automatically shuts off after fifteen minutes of continuous use at the same level. The Heat and Massage Option is powered by a low-voltage transformer with an eight-foot hospital-grade 120V cord. The Lumex Heat and Massage Unit meets UL60601-1 Standards for electric medical devices.



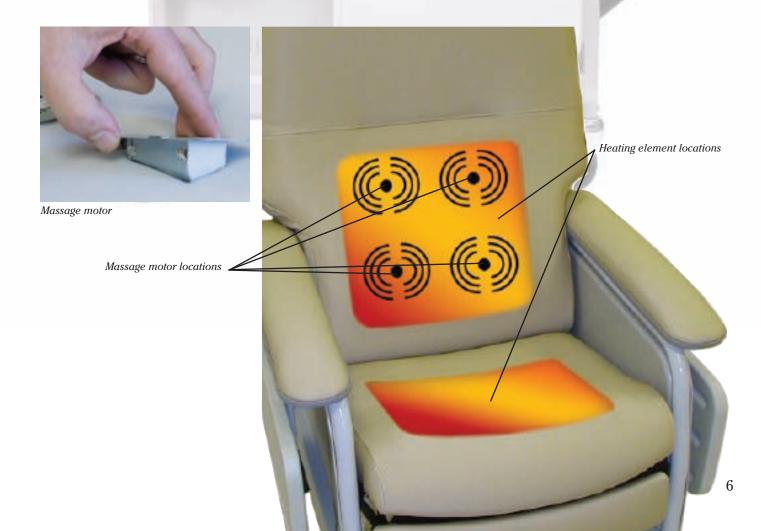



LUMEX FR566G CLINICAL CARE

The Lumex Deluxe Clinical Care Recliner sets a new standard of comfort for the patient and ease of use for the caregiver. The Lumex Deluxe Clinical Care Recliner is EXTRAordinary in many ways:

☐ The unique patent pending recline mechanism enables the patient to position the back from TV position to full recline. The wide, contoured seat and back, along with an integrated footrest, provide an enhanced level of comfort. Built-in arm pads in the back give

improved support for the arms while in Trendelenburg. Wall saver armrest caps increase durability of the armrests and protect walls from marring.

☐ Higher floor to back height enables the caregiver to raise the patient up with less effort, reducing the chance for a back injury.

☐ Extra deep Trendelenburg position steeper angle improves patient stabilization.

- ☐ Four Tente® swivel casters facilitate transport and maneuverability in confined areas. Left and right rear top-locking casters provide positive lock. Directional wheel lock on right rear caster reduces "crabbing" of recliner and enables the chair to be moved down hallways with less effort.
- ☐ Durable, high-impact side tables sit flush with chair arm to conveniently hold prep materials. Easy to actuate by caregiver.
- ☐ Extra heavy duty recliner mechanism.
- ☐ Weight bearing capacity: 350 Pounds.



Features LIQUICELL® Technology

Easy to achieve Trendelenburg position. No more cables or handles. Trendelenburg

position can be actuated from either side of the recliner by a convenient foot pedal.

Vinyl covered

side panels



Multiple Positions for Enhanced Patient Comfort



LUMEX FR566DG CLINICAL CARE

☐ Innovative, easy to actuate drop-arm feature. Self-storing, dual drop-arms are simple for the caregiver to activate and enables a safer, easier lateral transfer for the resident and caregiver.

☐ The unique recline mechanism enables the patient to position the back from TV position to full recline. The wide, contoured seat and back, along with an integrated footrest provide an enhanced level of comfort. Built-in arm pads in the back give improved support for the arms while in Trendelenburg. Wall saver armrest caps increase durability of the armrests and

protect walls from marring.

☐ Easy to achieve Trendelenburg position. No more cables or handles, Trendelenburg position can be actuated from either side of recliner by a convenient foot pedal. Higher floor to back height enables the caregiver to raise the patient up with less effort reducing the chance for a back injury.

☐ Extra deep Trendelenburg position—steeper angle improves patient stabilization.

☐ Easily removed "pop-out" side panels enable easy access to recliner for cleaning and maintenance.

- Four Tente® swivel casters facilitate transport and maneuverability in confined areas. Left and right rear top-locking casters provide positive lock. Directional wheel lock on right rear caster reduces "crabbing" of recliner and enables the chair to be moved down hallways with less effort.
- ☐ Durable, high-impact side tables sit flush with chair arm to conveniently hold prep materials. Easy to actuate by caregiver.
- ☐ Extra heavy duty recliner mechanism.
- □ Weight bearing capacity: 350 Pounds.

PERMA**BLOK**3® ADVANCED VINYL PROTECTION

Features **LIQUICELL®** Technology











Lumex FR577RG Clinical Care

- □ Welded steel frame is constructed of round tubing resulting in improved structural integrity of the frame.
- ☐ Improved wider and softer armrest for patient comfort.
- □ Softer seat and back package for improved patient comfort. Varied density seat back and redesigned seating position provides excellent comfort and support during long term sitting.
- ☐ Easy to activate gas cylinder mechanism for quick Trendelenburg activation and improved lift.
- Durable, high impact side tables sit flush with chair arm to conveniently hold prep materials. Easy to actuate by caregiver.
- ☐ Easily removed "pop-out" side panels enable easy access to recliner for cleaning and maintenance.
- ☐ Three-position, plus Trendelenberg.☐ Four Tente® swivel casters facilitate transport and maneuverability in confined areas. Left and right rear top-locking casters provide positive lock. Directional wheel lock on right rear caster reduces "crabbing" of recliner and enables the chair to be moved down hallways with less
- effort.

 ☐ Weight Bearing Capacity: 290 Pounds.



PERMABLOK^{3®} ADVANCED VINYL PROTECTION

Features LIQUIGELL® Technology

BACK PROFILE: FR577 Featuring specially-contoured high-density foam, the

new back profile enhances patient comfort and support.



Infinite Position Back Design





Infinite Position Back Design



LUMEX ELECTRIC BARIATRIC CLINICAL CARE RECLINER FR588W

The Lumex Electric Bariatric Recliner is the perfect solution for larger patients who require a wider chair. Composed of a heavy gauge steel frame, this chair has a maximum weight capacity of 700 pounds! A four-function hand pendant activates the motor, allowing caregivers to safely recline and raise patients without strain. Quick-release trigger mechanism on push handle allows chair to swiftly move into Trendelenburg position. ☐ Fully equipped with rechargeable battery pack for full operation away from a wall outlet. ☐ Four swivel casters facilitate transport and maneuverability in confined areas. ☐ Durable, high impact side tables sit flush with chair arm to conveniently hold prep materials. Easy to activate by caregiver. ☐ Legrest and backrest raise and lower independently, giving the operator complete control over positioning. ☐ Comes complete with attaching head pillow for extra support and added comfort. ☐ Limited lifetime warranty on frame. ☐ One year warranty on electronics. ☐ Three years warranty on all other components. □ Weight bearing capacity: 700 Pounds. Features LIQUICELL® Technology The AEGIS Microbe Shield used on

Lumex Long Term Care and Patient Room Seating works by chemically bonding itself to the treated surface and rupturing the membranes of any bacteria that it comes in contact with which inhibits the growth and spread of microbial growth.

Available in Four Designer Colors









Blue Ridge (427) Doe Skin (851)

Jade (857)

Rosewood (863)

CLINICAL CARE RECLINERS ACCESSORIES

Leg Extensions 5699G

- ☐ Fits FR566G, FR577RG, and FR587W
- \Box Fits all round tube frame designs.
- \square Adds 5" to height of recliner.
- ☐ Set of four extension legs.

IV Pole Attachment 5782RG

- ☐ Offers a tip-free means of administering enteral feedings or IV therapy.
- ☐ For use on most Lumex Ortho-Biotic Seating with Lumex Bed Socket IV Pole (Model 2101A) sold separately.





Universal Head Pillow FR5658

- ☐ Easy-on hook and loop attachment, the Universal Head Pillow can be easily positioned for individualized patient comfort.
- ☐ Fits on almost every model of recliner.

Vintage Gold FR5658401 Warm Taupe FR5658409 Blue Ridge FR5658427 Imperial Blue FR5658432 Royal Blue FR5658454 Doe Skin FR5658851 Jade FR5658857 Rosewood FR5658863

Refer to page 17 to view colors.



SHIPPING CARTONS

All Lumex Clinical Care Recliners are now shipped in specially-designed cartons to reduce shipping damage.
Casters come pre-installed.
Tool-free assembly of recliners is quick and easy.

Assembly instructions are included with every recliner.



Back-saver design makes unpacking easier by reducing lifting.







Folding Set Up Screens


Sets up in minutes to afford individual

privacy. Vinyl reinforced panels with heat-sealed hems on snap-out rods, lightweight aluminum hinged tubing frame. Use with casters (included) or easily convert to tipped leg bases (included). Base allows screen to extend to full width and adjusts for multiple configurations. Easily dismantled for storage. White panels are 27.5" x 56".

4297W Three-section screen 4298W Four-section screen

4299C Replacement casters, set of 4 **4299W** Replacement white panel



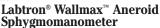


Large, 6" easy-to-read gauge is attached to a 110° swivel wall mounting bracket. Accumax™ cuff includes 8' coiled tubing. Cuff can be stored behind gauge.

222 Adult Blue

V222B Adult Latex-Free Black

Cuff and Faceplate
Replacement Crystal

Faceplate

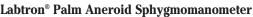

Durable, aluminum manometer with large, easy-to-read 65 mm dial. Integral latex-free, ribbed inflation bulb with dust filter. Deluxe push-button air release valve. Rotating bulb and valve provide ambidextrous use. Nylon cuff and latex-free one-tube bladder. Deluxe nylon carry case. Twenty year calibration warranty

215 Adult Grey
215C Child Grey
215X Large Adult Grey
216 Adult Orange



Labtron[®]Mobile Aneroid Sphygmomanometer

This mobile, height-adjustable aneroid sphygmomanometer features a large, 6" easy-to-read gauge mounted on a tip-resistant 5-wheel base. Accumax[™] cuff includes 8' coiled tubing. Cuff can be stored behind gauge.

223B Adult Black V223 Adult Black Latex-Free



Deluxe Revolving Stool With Foot Ring

Durable, 4" foam-padded, 14.5" diameter vinyl seat on ruggedly constructed tubular chrome-plated steel frame. Large base spread provides added stability. Footrest ring for comfort. Height adjustment 17.5" to 27.5".

1834 B Black seat

1834 C Casters for 1834B, set of 4



Heavy-Duty Footstool

- Rugged, welded, chrome-plated steel construction.
- Non-tip design.
- Platform topped with non-skid, ribbed mat.
- Equipped with reinforced rubber tips.
- 9.5"H x 15"W x 16"D
- Weight capacity 350 lb.

GF1841-2 2/cs

Safety Step-Up Stool

- Attractive chrome-plated finish.
- Non-tip design.
- Features 35"-high handle platform and steel-reinforced rubber tips.
- · Non-skid, ribbed rubber mat.
- 9.5"H x 15"W x 17"D
- Weight capacity 250 lb.

GF1840C-2 2/cs




Stainless Steel Stethoscope


The stainless steel dual head chestpiece has a non-chill ring for patient comfort. The internal spring binaural ensures acoustic sensitivity. 22" latex-free tubing keeps

acoustics intact, even in the noisiest settings. Comes complete with two pair of contoured

eartips. Lifetime warranty.


LAB-7100 Black, Adult
LAB-7100GY Grey, Adult
LAB-7200 Black, Pediatric
LAB-7300 Black, Infant
LAB-7300GY Grey, Infant

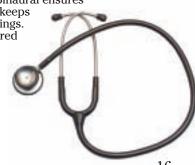

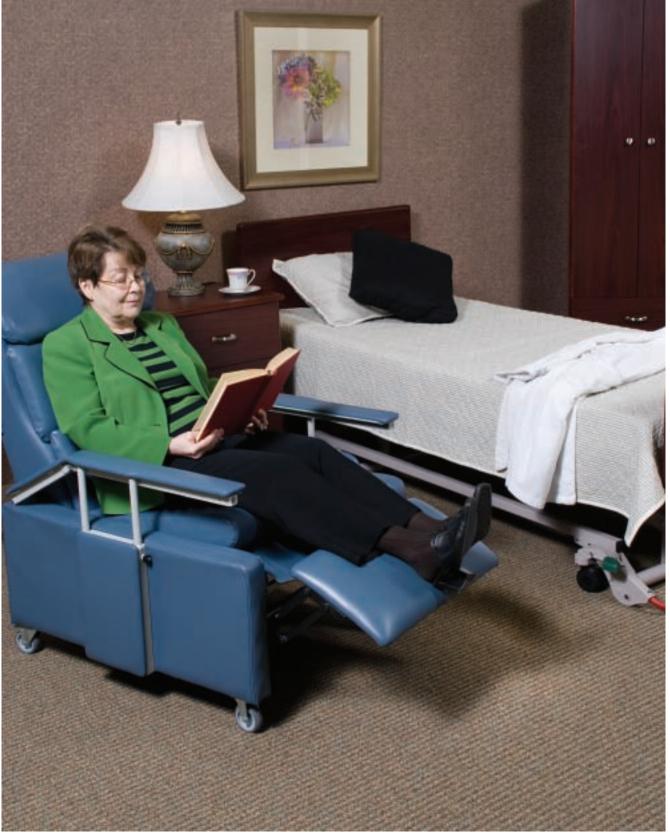



Vintage Gold (401) Warm Taupe (409) Blue Ridge (427) Imperial Blue (432)

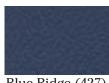









The AEGIS Microbe Shield used on Lumex Long TermCare and Patient Room Seating works by chemically bonding itself to the treated surface and rupturing the membranes of any bacteria that it comes in contact with which inhibits the growth and spread of microbial growth.

Jade (857)

Rosewood (863)

DROP-ARM RECLINERS

Lumex Drop-Arm Recliners FR3302

- ☐ Pillow back, three-positions, fire-rated upholstery.
- ☐ Combines the therapeutic benefits of position change with a more residential appearance. Designed to aid in circulation and to help ease discomfort from respiratory, heart and other ailments.
- ☐ Unique drop-arm mechanism facilitates lateral transfers. Easy to activate by caregiver.
- ☐ Ideal for knee and hip surgery patients.
- ☐ Three inch ball bearing casters ensure smooth transport; foot actuated swivel locking rear casters for a more
 - foot actuated swivel locking rear casters for a more secure lock.
- ☐ Unitized welded steel frame for durability. Flame and stain resistant vinyl upholstery available in a variety of colors.
- ☐ Weight Bearing Capacity: 250 Pounds.
- ☐ Four Fire Rated Upholstery Models are Available (Meets California Technical Bulletin 133 Flammability Standards)

Lumex Wrap-Around Warranty

- 1. Welded Frame: Limited Lifetime
- 2. Upholstery: Limited One Year
- 3. All Other Components: Limited Three Year

Available Colors:

Blue Ridge FR3302427 Doe Skin FR3302851 Jade FR3302857 Rosewood FR3302863

Refer to page 15 to view colors.

Lumex Drop Arm Recliner Accessories

- ☐ Anti-Tip Wheels 3306 One pair
- ☐ Tray Table 331000



Anti-Tip Wheels 3306



Tray Table 331000

THREE POSITION RECLINERS

Lumex Three Position Recliner 574G

- ☐ The Standard of the Industry and the recliner that started it all! This first recliner specifically designed to improve the comfort and mobility of residents in extended care facilities. The Lumex Three Position Recliner easily adjusts into any one of three positions: sitting, TV or full recline. Recliners offer numerous therapeutic benefits including position changes that can aid in circulation and ease discomfort from respiratory, heart and other ailments.
- ☐ A position lock allows only the caregiver to secure the chair in any one of the three easily attained positions. Recliner encourages natural relaxation and improved

☐ Flex-o-Lator seat support system increases comfort by reducing the risk of "bottoming out".

blood circulation.

☐ Ergonomically designed seat and back with varied density foam construction provides correct support. Articulating headrest quickly adjusts to maximize individual comfort requirements.

☐ For safer, easier transport, the Lumex Three Position Recliner features two four-inch swivel casters with wheel locks and four-inch fixed front wheels.

☐ The molded activity table can be fixed in position or left unsecured for easy exit from the recliner.

☐ Durable welded steel epoxy coated frame is easy to clean and maintain.

☐ Weight Bearing Capacity: 250 pounds.



Frame: Limited 2 Years
 Upholstery: Limited 1 Year

3. All Other Components: Limited 1 Year



Side Bolsters 599G

Available Colors:

Vintage Gold 574G401
Warm Taupe 574G409
Blue Ridge 574G427
Imperial Blue 574G432
Royal Blue 574G454
Doe Skin 574G851
Jade 574G857
Rosewood 574G863

Fire Rated Colors

Meets California Technical Bulletin 133 Flammability Standards

Vintage Gold Warm Taupe Blue Ridge Imperial Blue Royal Blue FR574G432 FR574G432 FR574G454 Doe Skin Jade FR574G851 FR574G857 Rosewood FR574G863

20

Preferred Care® Recliners Series

The Lumex Preferred Care Recliner Series offers facilities the widest selection of sizes and features of any specialty healthcare seating. When it comes to quality care, one size does not fit all and proper fit is an important consideration. Lumex Preferred Care Recliners are available in regular and extra-wide versions. Based upon the original 565G, the 565WG offers the same features scaled to accommodate the resident's unique proportions.

Standard Features on all Lumex Preferred Care Recliners:

- ☐ Infinite recliner back patented mechanism provides infinite positioning of the back, ranging from 95 degrees to 145 degrees for increased resident comfort and pressure relief. The reclining back can be locked in any position by the caregiver.
- ☐ A unique suspension system self-adjusts to each person, cradling the resident from head to heel and virtually eliminating the risk of "bottoming out".
- ☐ Extra-wide padded armrests and side panels maximize comfort.
- ☐ Durable epoxy coated welded steel frame is easy to clean and maintain.
- ☐ Quick-release upholstery allows for easy cleaning, removal and replacements. Vinyl upholstered back is sulfide stain, mildew and oil resistant.

565G

- ☐ A variety of accessories are available to customize the chair for specific individual needs.
- ☐ Lumex Wrap-Around Warranty
 - 1. Welded Frame: Limited Lifetime
 - 2. Upholstery: Limited One Year
 - 3. All Other Components: Limited Three Year
- ☐ Weight Bearing Capacity:

565DG

- 1. 565G, 565TG, and 565DG: 250 lb.
- 2. 565WG: 400 lb.

Preferred Care® Recliners Series

Lumex[®] **Preferred Care**[®] **Recliner 565G**

The Lumex 565 Patented recline mechanism allows for infinite positioning of recline. Reduced armrest and seat height enables the recliner to fit under dining/activity tables and allows resident to "scoot" and more easily reposition chair by themselves. Two swivel casters with two rear locking casters enable the Preferred Care Recliner to be easily maneuvered in small spaces and moved down hallways.

Available Colors:

Vintage Gold 565G401 Warm Taupe 565G409 Blue Ridge 565G427 Imperial Blue 565G432 Doe Skin 565G851 Jade 565G857 Rosewood 565G863

Fire Rated Colors

Meets California Technical Bulletin 133 Flammability Standards

Blue Ridge FR565G427
Jade FR565G857
Rosewood FR565G863
Refer to page 17 to view colors.



Shown With Optional

Accessories.

Lumex[®] Extra-Wide Preferred Care[®] Recliner 565WG

All of the features and benefits of the Lumex Preferred Care Recliner in an extra-wide width. The 24" width is ideal for residents who require a wider seat width. Directional wheel lock on right rear caster reduces "crabbing" of recliner and enables the chair to be moved down hallways with less effort. Weight bearing capacity: 400 lb.

Available Colors:

Vintage Gold 565WG401 Warm Taupe 565WG409 Blue Ridge 565WG427 Imperial Blue 565WG432 Doe Skin 565WG851 Jade 565WG857 Rosewood 565WG863

Fire Rated Colors

Meets California Technical Bulletin 133 Flammability Standards

Blue Ridge FR565WG427
Jade FR565WG857
Rosewood FR565WG863
Refer to page 17 to view colors.





PREFERRED CARE® RECLINER ACCESSORIES

Activity Tray Table 5644G & 5644WG

- Strong polypropylene tray table offers a supportive, stable work surface for activities.
- ☐ Includes receiving arms, storage bracket, and strap.
- ☐ Easily stored in the back of chair.
- ☐ Patent pending "lock/ lock-out" key gives flexibility to the caregiver.
- 1 5644G fits 565G, 565DG and 565TG
- □ 5644WG fits 565WG



Foot-Drop Assembly 5657G

- ☐ Helps keep resident's feet positioned correctly to help avoid foot-drop.
- □ Note: Foot-Drop Assembly will not fit 565WG.

Colors Available:

Blue Ridge 5657G427 Jade 5657G857 Rosewood 5657G863

Refer to page 17 to view colors.



Head and Body Bolster 56461, 56462, and 56471

- □ Soft, comfortable positioning devices.
- ☐ Available separately, they are easily installed without tools.
- ☐ Helps to maintain proper body alignment.
- □ Note: Body Bolster will not fit 565WG.

Head Bolster (56461)

Colors Available:

| 9 |
|---|
| 2 |
| 1 |
| 1 |
| |

Fire Rated Colors

Meets California Technical Bulletin 133 Flammability Standards

| Vintage Gold | FR56461401 |
|---------------|------------|
| Warm Taupe | FR56461409 |
| Blue Ridge | FR56461427 |
| Imperial Blue | FR56461432 |
| Royal Blue | FR56461454 |
| Doe Skin | FR56461851 |
| Jade | FR56461857 |
| Rosewood | FR56461863 |

Body Bolster (56471)

Colors Available:

| Vintage Gold | 56471401 |
|--------------|----------|
| Warm Taupe | 56471409 |
| Royal Blue | 56471454 |
| Doe Skin | 56471851 |
| Rosewood | 56471863 |

Fire Rated Colors

Meets California Technical Bulletin 133 Flammability Standards

| Vintage Gold | FR56471401 |
|---------------|------------|
| Warm Taupe | FR56471409 |
| Blue Ridge | FR56471427 |
| Imperial Blue | FR56471432 |
| Royal Blue | FR56471454 |
| Doe Skin | FR56471851 |
| Jade | FR56471857 |
| Rosewood | FR56471863 |

Refer to page 17 to view colors.

Head Bolster for 565WG (56462)

Colors Available:

| Vintage Gold | 56462401 |
|---------------|----------|
| Warm Taupe | 56462409 |
| Blue Ridge | 56462427 |
| Imperial Blue | 56462432 |
| Royal Blue | 56462454 |
| Jade | 56462857 |

Fire Rated Colors

(Meets California Technical Bulletin 133 Flammability Standards)

| Vintage Gold | FR56462401 |
|---------------|------------|
| Warm Taupe | FR56462409 |
| Blue Ridge | FR56462427 |
| Imperial Blue | FR56462432 |
| Royal Blue | FR56462454 |
| Doe Skin | FR56462851 |
| Jade | FR56462857 |
| Rosewood | FR56462863 |
| | |



Preferred Care® Recliner Accessories



Adjustable Footrest LX5640G

- ☐ Enhances patient positioning, comfort and aids in the prevention of foot-drop.
- Footrests tip up and out of way to allow access for easy patient transfer.
- ☐ Attaches without the use of tools.
- ☐ For use on Preferred Care Recliner.

Universal IV Pole Holder 5782RG

- ☐ Offers a tip-free means of administering internal feedings or IV therapy.
- For use on most Lumex Ortho-Biotic Seating with Lumex Bed Socket IV Pole (2101A) sold separately.



Drainage Bag Hook with Privacy Pouch 5645

- ☐ Helps improve resident dignity by disguising the drainage bag in a pouch that matches the upholstery.
- ☐ For use on 565 Preferred Care Recliners only.

Available Colors:

Blue Ridge 5645427 Jade 5645857 Rosewood 5645863

Refer to page 17 to view colors.



HIP CHAIR GF4404

The Lumex Everyday Hip Chair is ideal for people who have arthritis or for patients who have undergone hip or knee surgery. The 27" high seat helps to reduce stress and fatigue on joints when entering and exiting the chair. Gently curving arms on the frame make sitting and standing easier. Padded armrests and a footrest provide greater comfort while sitting. The anti-slip strip on the footrest provides added safety. Sturdy metal frame.

Weight Bearing Capacity: 350 Pounds

Warranty:

Frame: Two Years Upholstery: One Year

Colors Available:

Blue Ridge GF4404427 Jade GF4404857 Rosewood GF4404863

Refer to page 17 to view colors.



LUMEX SPECIALTY HEALTHCARE SEATING: CUSTOMER SUPPORT AND SERVICE

All Lumex Specialty Healthcare Seating is protected by our **Wrap-around Limited Warranty.**

Lumex Specialty Healthcare Seating purchases are not just about price. They are about COST IN USE.


Cost in Use is a concept that speaks to the Life Cycle of the product, not just the initial capital expense.


Lumex Specialty Healthcare Seating is your solution because service parts and support are available to help you reduce the cost in use of your healthcare seating. All upholstery components can easily be replaced. This greatly extends the life of your Lumex Seating Product.

The Graham-Field Website, www.grahamfield.com, is always available and is your resource for:

- ☐ Lumex Specialty Seating Schematic Drawings that clearly detail replacement parts and components. Schematics can be found on the website under the product Additional Information tab.
- ☐ Live Chat with a Graham-Field Technical Support Specialist
- ☐ A wide array of Specialty Seating accessories



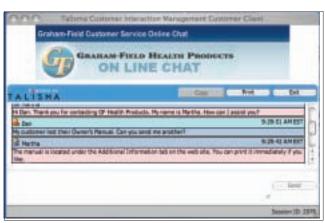

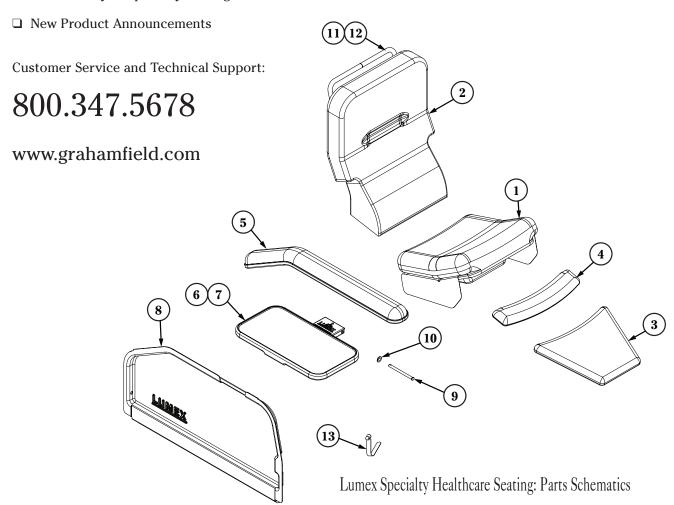

Specifications and Dimensions

| GF4404 42" | 565DG 48" | 565TG 48" | 565WG 48" | 565G 48" | 574G 47" | FR3302 44" | FR588W 55" | FR587W 50" | FR577RG 47" | FR566DG 51 ½" | FR566G 50" | Item Overall Height |
|-------------------|------------------|-------------------|------------------|-----------------|-----------------|-------------------|-------------------|------------|-------------|---------------|-------------------|---|
| 2" 25 ¾" | 3" 28" | 3" 28" | 32" | 8" 28" | 7" 25" | 4" 31" | 5" 38" |)" 33" | 7" 31" | 1/2" 33" |)" 31" | rall Overall ght Width |
| 1" 26" | 36" | 36" | 37" | 36" | 41" | 38" | 59" | 36" | 35" | 38" | 43" | all Overall h Depth |
| 18" | 19 ½" | 19 ½" | 24" | 19" | 20" | 23" | 25 1/4" | 23" | 22" | 23" | 21" | Width Between Arms |
| 31" | 23 ½" - 29 ½" | 23 1/2" - 29 1/2" | 28" | 28" | 26" | 19 1/2" - 25 1/2" | 32" | 26 ½" | 25 ½" | 28" | 27 ½" | Arm Height From Floor |
| 6" | 0" - 8 ½" | 0" - 8 1/2" | 8 1/2" | 8 ½" | 9" | 5" | 8" | 8" | 8" | 8" | 8" | Arm Height From Seat |
| 26" | 22" | 22" | 20" | 20" | 21" | 19" | 25 ½" | 21 ½" | 22" | 22" | 22" | Seat Height |
| 18" | 19" | 19" | 24" | 19" | 19" | 21 ½" | 25 ¾" | 23 ½" | 21 ½" | 20 ½" | 20 ½" | Seat Width |
| 18" | 19 ½" | 19 1/2" | 19 ½" | 19 ½" | 19" | 20" | 24" | 20 1/2" | 21" | 19" | 20" | Seat Depth |
| _ | | _ | 1 | | _ | _ | 71" | 76" | 74" | 76" | 75" | Depth Fully Reclined In Trendelenburg |
| _ | | _ | 1 | | _ | - | 58" | 53" | 51" | 53" | 51" | Width With Both Side Tables in Upright Position |
| 52 lb | 85 lb | 90 lb | 91 lb | 80 lb | 110 lb | 109 lb | 308 lb | 138 ½ lb | 127 lb | 135 lb | 145 lb | Shipping Weight |
| 350 lb | 250 lb | 250 1ь | 400 lb | 250 lb | 250 lb | 250 lb | 700 lb | 450 lb | 290 lb | 350 lb | 350 lb | Max. Weight Capacity |

For the most up-to-date information and product specifications, please visit our website www.grahamfield.com.



GF Health Products, Inc. Atlanta, Georgia 30360 800.347.5678 www.grahamfield.com